

Hemostatic Granules(Premium/Standard)

Stop Severe Bleeding
Mould to the complex shaped wounds

Premium Grade

- 1.No heat - no burning
- 2.Strong wound adhesion and easy debridement
- 3.Rapid Clotting: Fluid absorption within 30s ≥ 10 times
- 4.Similar to human physiological pH, no stimulation
- 5.Works in hypothermic conditions and anti-coagulant conditions(such as Coumadin and Heparin)
- 6.Bio-material , easy degraded, the remaining will be safely absorbed

Standard Grade

- 1.No heat - no burning
- 2.Strong wound adhesion
- 3.Rapid Clotting: Fluid absorption within 60s ≥ 10 times
- 4.Works in hypothermic conditions and anti-coagulant conditions(such as Coumadin and Heparin)




arterial bleeding / venous bleeding /
penetrating wounds /
superficial wounds / deep wounds
complex shaped wounds
open fractures accompanied by bleeding



When severe bleeding, every second counts
think ChitoGuard™


How to use Chitosan Hemostatic Granules

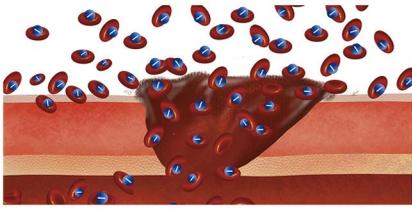

1. Open Wound




1. Identify the wound site and clean it quickly.

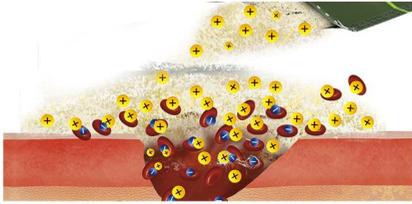

2. Use ChitoGuard™ Granules


2. Release enough granules and cover with gauze.

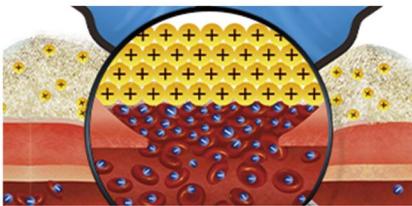

3. Chitosan coagulates blood



3. Apply firm pressure directly to the wound until bleeding stops.

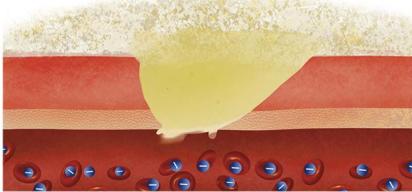

4. Clotted

4. Bandage and seek urgent medical attention.



ChitoGuard™ Products

Standard Grade (Chitosan)	Granules	1g/2g/3g/6g/9g/15g/20g/25g...
	Granules Applicator	6g
	Granules Xquick Plug	5g/6g
	Gauze Pad	7.5*10cm/7.5*7.5cm/8*8cm/10*10cm...
	Z-fold Gauze	7.5*150cm/7.5*180cm/7.5*300cm...
Premium Grade (CM-Chitosan)	Rolled Gauze	7.5*150cm/7.5*180cm/7.5*300cm...
	Granules	1g/2g/3g/6g/9g/15g/20g/25g...
	Granules Applicator	6g
	Granules Xquick Plug	5g/6g

Qualification



Contact us:

Telephone: (+86)-532-58203680

Email: sales@chitoguard.com

Website: www.chitoguard.com

Address: No.10, Shengshui Rd, Laoshan District, Qingdao, China

EN ISO13485:2016

**When severe bleeding, every second counts
think ChitoGuard™**